

PRODUCT INFORMATION

Tag	C-Flag&Strep Tag
Target	CHRM3
Synonyms	EGBRS; HM3; PBS
Description	Human CHRM3-Strep full length protein-synthetic nanodisc
Delivery	In Stock
Uniprot ID	P20309
Expression Host	HEK293
Protein Families	Druggable Genome, GPCR, Transmembrane
Protein Pathways	Calcium signaling pathway, Neuroactive ligand-receptor interaction, Regulation of actin cytoskeleton
Molecular Weight	The human full length CHRM3-Strep protein has a MW of 66.1 kDa
Formulation & Reconstitution	Lyophilized from nanodisc solubilization buffer (20 mM Tris-HCl, 150 mM NaCl, pH 8.0). Normally 5% - 8% trehalose is added as protectants before lyophilization. Please see Certificate of Analysis for specific instructions. Do not use solvents with a pH below 6.5 or those containing high concentrations of divalent metal ions (greater than 5 mM) in subsequent experiments.
Storage & Shipping	Store at -20°C to -80°C for 12 months in lyophilized form. After reconstitution, if not intended for use within a month, aliquot and store at -80°C (Avoid repeated freezing and thawing). Lyophilized proteins are shipped at ambient temperature.
Background	The muscarinic cholinergic receptors belong to a larger family of G protein-coupled receptors. The functional diversity of these receptors is defined by the binding of acetylcholine and includes cellular responses such as adenylate cyclase inhibition, phosphoinositide degeneration, and potassium channel mediation. Muscarinic receptors influence many effects of acetylcholine in the central and peripheral nervous system. The muscarinic cholinergic receptor 3 controls smooth muscle contraction and its stimulation causes secretion of glandular tissue. Alternative promoter use and alternative splicing results in multiple transcript variants that have different tissue specificities.
Usage	Research use only
Conjugate	Unconjugated



ELISA assay to evaluate CHRM3-Strep-Nanodisc
0.2µg Human CHRM3-Strep-Nanodisc per well

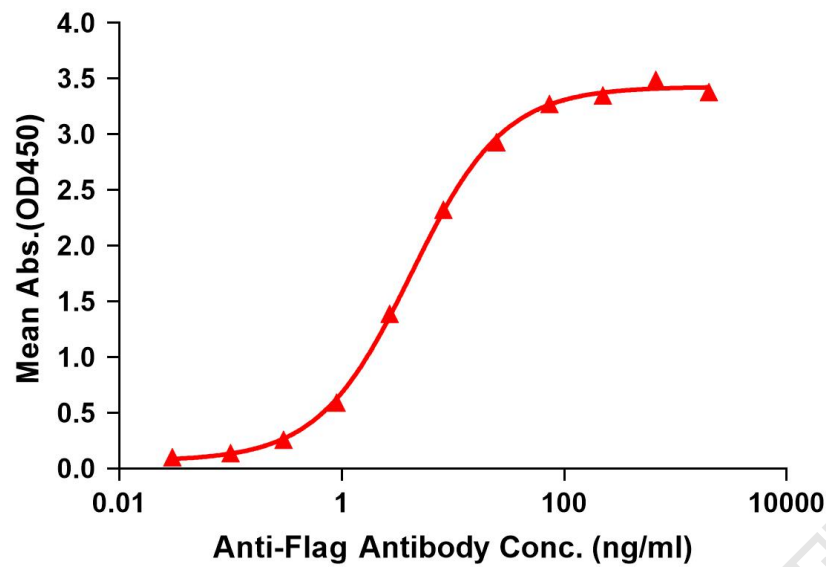


Figure 1. Elisa plates were pre-coated with C-Flag&Strep Tag CHRM3-Strep-Nanodisc (0.2µg/per well). Serial diluted anti-Flag monoclonal antibody solutions were added, washed, and incubated with secondary antibody before Elisa reading. From above data, the EC50 for anti-Flag monoclonal antibody binding with CHRM3-Strep-nanodisc is 4.223ng/ml.



Figure 2. Human CHRM3-Strep-Nanodisc, C-Flag&Strep Tag on SDS-PAGE

